

FRED EASYPORT PLUS

Small and powerful defibrillator - for anyone, anytime, anywhere




FRED EASYPORT PLUS

The ideal companion on the go

SCHILLER's pocket-size defibrillator FRED easyport plus:

- 780 grams light, 46x150x143 mm (hxlxw)
- Available as AED First, Auto or Manual
- CPR feedback options

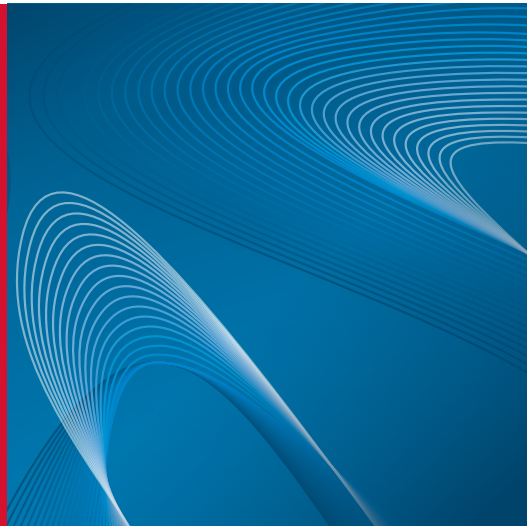

EASY TO CARRY

As a pocket-sized defibrillator, the FRED easyport plus weighs only 780 grams and combines small size with great functionalities.

It comes with a standard or a small carrying bag, depending on which accessories need to be packed with the defibrillator.

SHOCK AS SOON AS IT IS REQUIRED

FRED easyport plus features a quick analysis system and charges the required energy systematically so that the shock can be delivered as soon as the device advises it.

CPR FEEDBACK

Defibrillation alone is not enough. To maintain a good supply of oxygen to the organs, cardiopulmonary resuscitation (CPR) must be performed correctly to restore proper blood flow.

FRED easyport plus helps you achieve highest efficiency in chest compressions by:

- Indicating the correct depth, rate and recoil
- Setting a regular rhythm through the metronome

When using the ARGUS LifePoint sensor, the rescuer receives audio and visual assistance via colour-coded on-screen messages during chest compressions.

FLEXIBLE BATTERY CONCEPT

It can be operated with a rechargeable or a disposable battery. The external charger can charge two batteries simultaneously.



INFORMATION TRACKING

FRED easyport plus records intervention data, allowing doctors and emergency physicians to gain a clearer overview of the patient's condition.

- ❖ CPR quality
- ❖ ECG curve
- ❖ Impedance curve
- ❖ Events

INFORMATION RETRIEVAL

Thanks to Wi-Fi and USB functionalities, saved files can be easily reviewed and managed by SCHILLER's data management solution SEMA.

FIND THE DEVICE THAT SUITS YOU

FRED easyport plus is available in three versions to fit the different user levels.



AED FIRST

FRED easyport plus semi-automatic AED with its grey cover requires the rescuer's confirmation to deliver a shock.

An ECG curve can be displayed on the large LCD colour screen.



AED AUTO

The fully automatic FRED easyport plus requires no action from the rescuer. After the heart rate analysis, the device determines whether a shock is necessary. If it is, the shock is delivered automatically.

An ECG curve can be displayed on the large LCD colour screen.



AED MANUAL

Distinguishable by its red cover, the manual FRED easyport plus allows doctors and paramedics to switch to manual override. Thanks to the ECG display, healthcare professionals can determine if a shock is necessary.



SCHILLER
FRED easyport plus

00:00:04
PLUG AND APPLY ELECTRODES

ELECTRODES
FRED PA-1
RFID
ADULT



SMALL AND POWERFUL

FRED easyport plus weighs only 780 grams and delivers shocks of up to 200 joules. The Multipulse Biowave technology ensures an effective and safe defibrillation shock while minimising damage to the myocardium.

GAIN PRECIOUS SECONDS

With pre-connected defibrillation electrodes, heart rate analysis starts right away to reduce time until the shock is delivered.

ECG SIGNAL MODE

The 2-lead ECG signal mode is used when defibrillation is not indicated but the patient's cardiac rhythm is to be observed by medical staff.

ALWAYS READY

FRED easyport plus performs regular automatic self-tests to ensure it works correctly. Its status is indicated by the flashing on/off button.

DATA TRANSMISSION TO SEMA

Intervention data can be sent via Wi-Fi or USB for review by healthcare professionals.

DESIGNED TO LAST

With IP44 classification, FRED easyport plus resists splashing water.





POCKET-SIZED DEFIBRILLATOR

FRED easyport plus measures only 150 x 143 millimetres – this photo shows its actual size (1:1)!

SUITABLE FOR ALL ENVIRONMENTS

FRED easyport plus has a large colour LCD screen that automatically adjusts its brightness to help visualise indications at a glance.

INTUITIVE PAEDIATRIC MODE

While using the same electrodes, just press this button to change the patient mode.

CPR FEEDBACK

The integrated metronome sets the correct rhythm and guides the rescuer during the intervention.

In addition, the ARGUS LifePoint sensor assists the rescuer with audible and visual cues on screen to ensure effective chest compressions.

OPTIMISED USER INTERFACE

Three simple steps with illustrations and spoken instructions for a clear understanding.



Americas
SCHILLER Americas Inc.
Doral, Florida 33172

North America:
Phone +1 786 845 06 20
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com

Latin America & Caribbean
Phone +1 305 591 11 21
Fax +1 786 845 06 02
sales@schilleramericas.com
www.schilleramericas.com



France
SCHILLER Médical
F-67160 Wissembourg
Phone +33 3 88 63 36 00
info@schiller.fr
www.schiller-medical.fr



France (distribution France)
SCHILLER France S.A.S.
F-77608 Bussy St Georges
Phone +33 1 64 66 50 00
contact@schillerfrance.fr
www.schiller-france.com



Serbia
SCHILLER d.o.o.
11010 Beograd
Phone +381 11 39 79 508
info@schiller.rs
www.schiller.rs



Slovenia
SCHILLER d.o.o.
2310 Slovenska Bistrica
Phone +386 2 843 00 56
info@schiller.si
www.schiller.si



Asia
SCHILLER Asia-Pacific / Malaysia
52200 Kuala Lumpur, Malaysia
Phone +603 6272 3033
sales@schiller.com.my
www.schiller-asia.com



Germany
SCHILLER Medizintechnik GmbH
D-85622 Feldkirchen b. München
Phone +49 89 62 99 81 0
info@schillermed.de
www.schillermed.de



Spain
SCHILLER ESPAÑA, S.A.
ES-28232 Las Rozas/Madrid
Phone +34 91 713 01 76
schiller@schiller.es
www.schiller.es



Austria
SCHILLER Handelsgesellschaft m.b.H.
A-4040 Linz
Phone +43 732 70 99 0
sales@schiller.at
www.schiller.at



India
SCHILLER Healthcare India Pvt. Ltd.
Mumbai – 400 059, India
Phone +91 22 6152 3333 / 2920 9141
sales@schillerindia.com
www.schillerindia.com



Switzerland
SCHILLER Schweiz AG
CH-8912 Obfelden
Phone +41 44 744 30 00
sales@schiller-schweiz.ch
www.schiller-schweiz.ch



China
Alfred Schiller (Beijing) Medical Technology
Co., Ltd.
100102 Beijing, China
Phone +86 10 84565453
info@schillerchina.com
www.schillermedical.cn



Poland
SCHILLER Poland Sp. z o.o.
PL-02-729 Warszawa
Phone +48 22 843 20 89 / +48 22 647 35 90
schiller@schiller.pl
www.schiller.pl



Turkey
SCHILLER TÜRKIYE
Okmeydani-Sisli – Istanbul
Phone +90 212 210 8681
info@schiller.com.tr
www.schiller-turkiye.com



Croatia
SCHILLER d.o.o.
10000 Zagreb
Phone +385 1 309 66 59
info@schiller.hr
www.schiller.hr



Russia & C.I.S.
AO SCHILLER.RU
119049 Moscow, Russia
Phone +7 (495) 970 11 33
mail@schiller.ru
www.schiller.ru



United Kingdom
SCHILLER UK Ltd.
Bromborough, Wirral,
Merseyside, CH62 3RJ
Phone +44 1333 312150
sales@schilleruk.com
service@schilleruk.com
www.schilleruk.com

SCHILLER reserves the right to make changes in design, specifications and features shown or discontinue the product described at any time without notice or obligation. Images may be representative. Some products may not be regulatory cleared or released for sale in all markets. Contact your SCHILLER representative for the most current information and availability.

